

ANZEN™ safety scalpel

Instructions For Use

1 Attaching Blade Cartridge

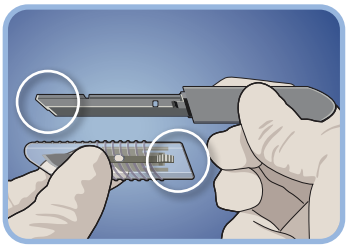

Align Handle and Cartridge orientations by matching end angles.

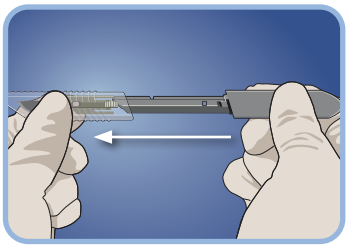


Insert Handle into Cartridge.

2 Extending Blade




Depress Button on Cartridge.

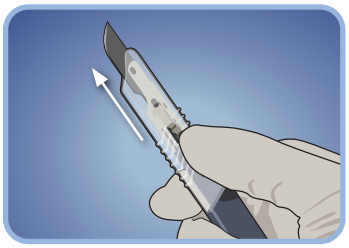


Slide Button forward until it locks.

3 Retracting Blade

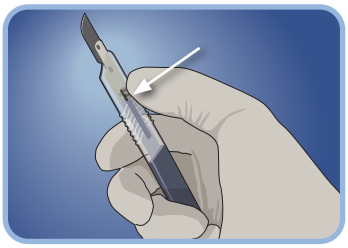


Depress Button on Cartridge.

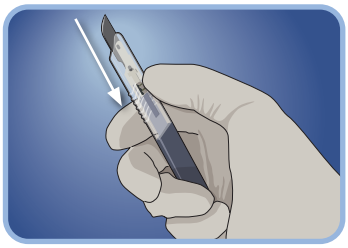


Slide Button back until it locks.

4 Detaching Blade Cartridge

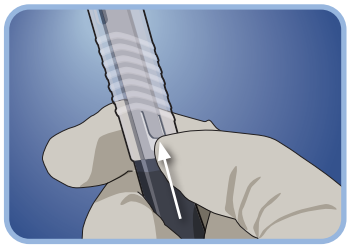


Lift Release Tab (Dorsal View).

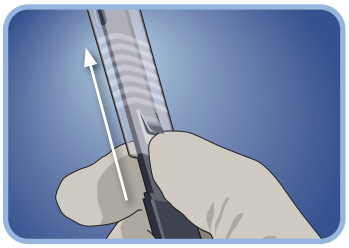


Slide cartridge forward to eject.

Cartridge is single-use only. Reuse may result in patient cross-contamination. After use, discard into a sharps container.

Ventral View (Button Side)

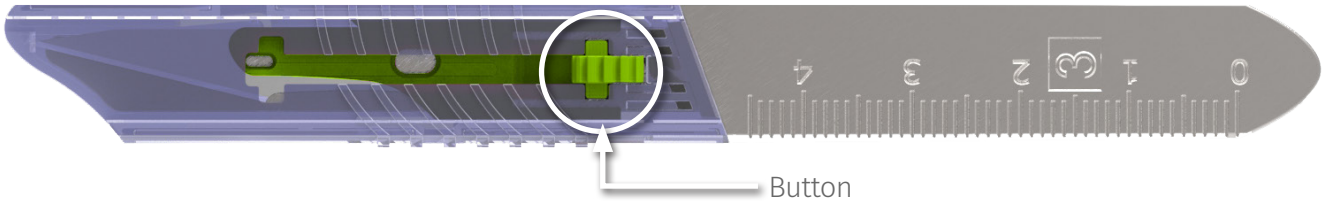

Dorsal View (Release Tab Side)

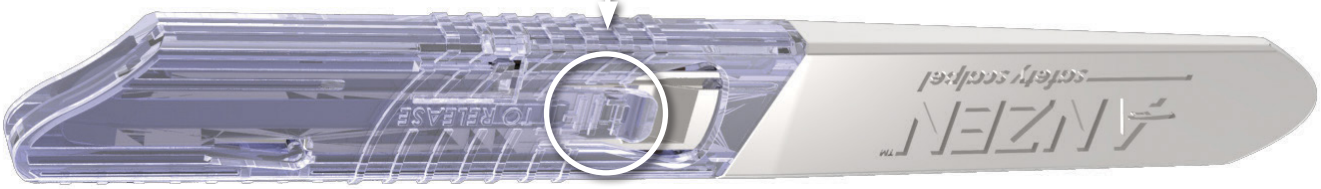

1. **Warnings:** For single patient use only. Reuse may result in patient cross-contamination.
2. Do not use if the package is broken
3. Do not store in high temperature, high humidity environment or with corrosive gas.
4. Do not use if product is damp or expired.
5. The **Anzen™ Safety Scalpel** is designed to be used by medical professionals.
6. After use, discard into a sharps container.
7. The **Anzen™ Safety Scalpel** handle can be autoclaved and reused as indicated in the instructions for use available using the website link.
8. In the event of a serious incident, contact MediPurpose at sales@medipurpose.com and/or the competent authority of the Member State in which the incident occurred.

Visit anzenscalpel.com or contact us:

MEDIPURPOSE PTE. LTD.
10 Anson Road, #12-08
International Plaza, Singapore 079903

EU REP Advena Ltd.
Tower Business Centre, 2nd Flr.,
Tower Street, Swatar, BKR 4013 Malta

CE 0197 MD

2862 Buford Hwy, Suite 105, Duluth, Georgia 30096, USA
tel +1-404-855-3655 | fax +1-770-559-4382 | sales@medipurpose.com

A product of MediPurpose™
www.medipurpose.com

©2026 MediPurpose Pte Ltd